



Primary Cell KO Pools:

Advance your immunology research with edited-to-order primary cells

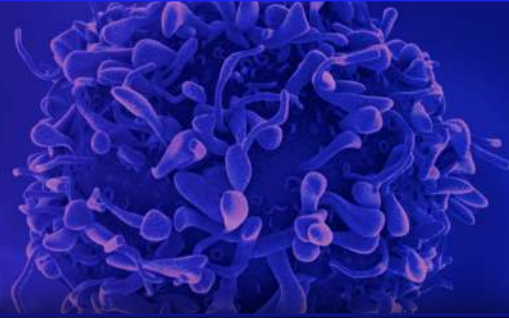

High editing efficiency with reliable post-editing functionality. Order Primary Cells supplied by EditCo or work with our team to edit your primary cell type of interest.

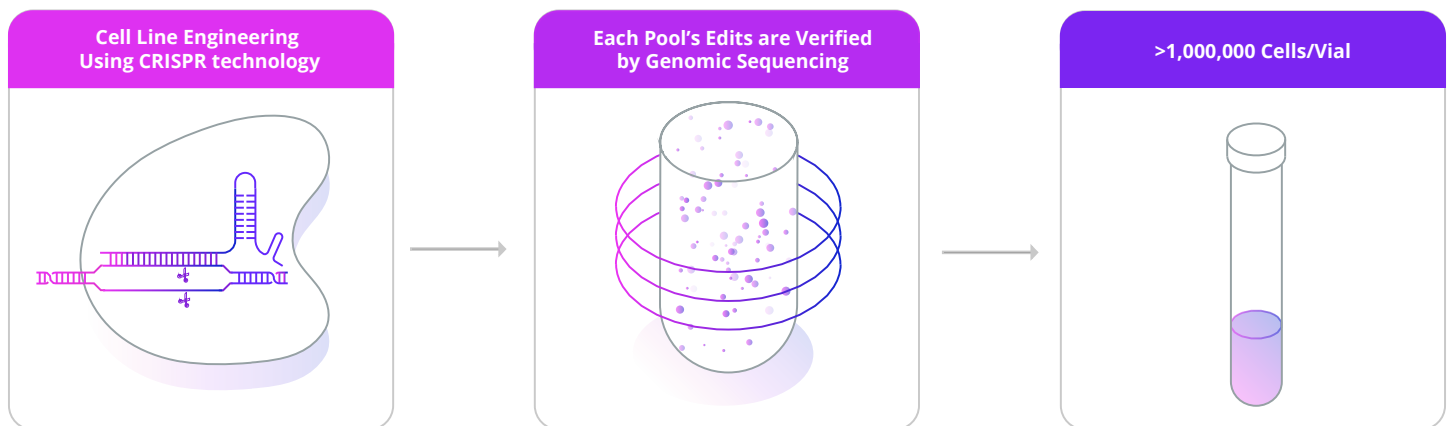

OVERVIEW

Skip to the main experiments with edited T cells

- ✓ **Confident Results:** 80% editing efficiency guaranteed smart multi-guide design.
- ✓ **Faster Results:** A 7-day editing protocol delivers results in 2 weeks or faster.
- ✓ **Flexibility:** Choose from EditCo-supplied cells or onboard your T cells.
- ✓ **Functionality:** High cell viability and editing levels without affecting functionality.

T cells: Key Product Deliverables

- ✓ Regular updates on your order's progress
- ✓ 2 vials of edited cell pools with >1,000,000 cells/vial (pools consist of a heterogeneous population of edited and unedited cells)
- ✓ Control-transfected cell pools (2 vials)
- ✓ Sequence of synthetic gRNA used
- ✓ Primer sequences used for PCR and Sanger Sequencing
- ✓ ICE Sanger sequencing analysis reports for each edited clone after expansion
- ✓ Comprehensive QC report that includes the following information: mycoplasma test (positive/negative) and passage number

Using other primary cell lines?
Reach out to us about custom projects!

Scan to learn more




Efficient genome editing: EditCo guarantees no less than 80% editing efficiency and >70% cell viability at time of freeze for both edited and CD8+ T cells.

CD4+ T Cells

Donor ID	Donor 1	Donor 2	Donor 3
Race	White	White	White
Smoker	Non-Smoker	Non-Smoker	Non-Smoker
Age	59	24	52
Gender	F	M	F

Table 1. Characteristics of edited T cell donors. Donors are screened for all major pathogens such as HIV, CMV and others.

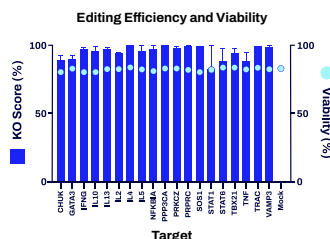


CD4+ T Cells Viability and Editing efficiency. Values represent the average of 3 donors. Error bars are S.E.M.

CD8+ T Cells


Donor ID	Donor A	Donor B	Donor C
Race	White	AA	White
Smoker	Non-Smoker	Non-Smoker	Non-Smoker
Age	29	29	23
Gender	F	F	M

Table 2. Characteristics of CD8+ T cell donors. Donors are screened for all major pathogens such as HIV, CMV and others.

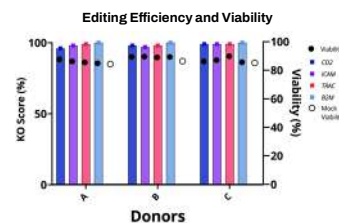

CD8+ T Cells Viability and Editing efficiency. Viability and editing efficiencies across all KO edits were measured in CD8+ T cell pools.

Cell fitness is not affected by the editing process. EditCo T cells maintain high cell viability and editing levels after thawing and expansion.

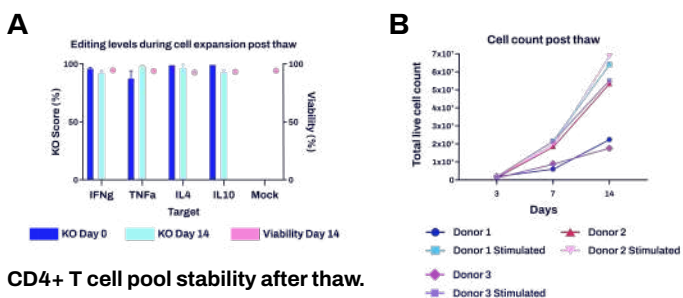


CD4+ T cell pool stability after thaw.


A. Pools of 4 different edits from 3 separate donors showed >90% viability with no decrease in editing efficiency after 14 days in culture.

B. Total live cell count over a 14 day period from 3 separate donors expanded with CD3/CD28 stimulation (stimulated) or IL2 alone.

C. PD1 expression in TCR activator-treated cultures peaked on day 3 in all tested donors (edited and mock samples). Isotype control staining over the same period is shown for reference.

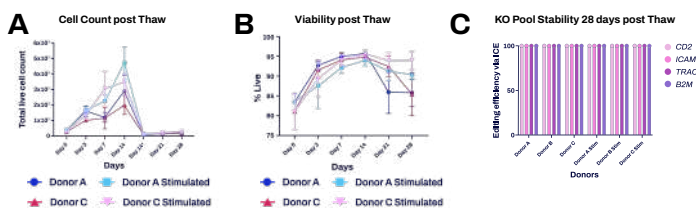


CD8+ T cell pool stability after thaw. **A.** Cell numbers and **B.** Viability were measured at the indicated time points. (*) Cell number was reduced to 10^6 /well on Day 14 to avoid overcrowding. **C.** Editing stability was measured at Day 28 in all cultures by PCR/Sanger sequencing. Charts show average values and Standard Error.

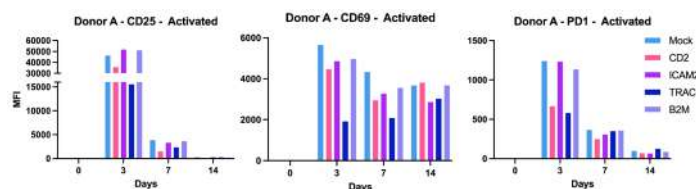


Transient expression of activation and exhaustion markers. Median Fluorescence Intensity (MFI) of activation markers (CD25 and CD69) or exhaustion markers (PD-1) in cultures of edited CD8+ T cells after thaw. Charts show average MFI and Standard Error.

Functional cells following editing. Edited CD4+ T cells produced levels of IFNg, IL2 and TNFa after mitogenic stimulation correlating with known target gene function. CD8+ T cells demonstrated antigen specific CD107a mobilization following CEF peptide stimulation.

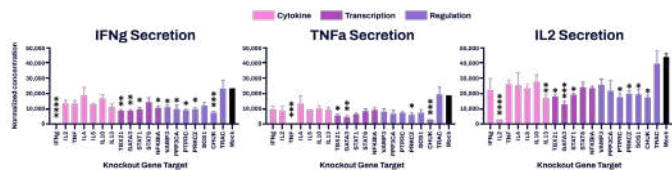

Cytokine secretion in edited cell pools following CD4+ T cell expansion. Edited cells from 3 donors listed in Table 1 were thawed and then stimulated with PMA/ionomycin for 6 hours. Supernatants were harvested and measured for concentrations of IFNG (left panel), IL2 (center panel) and TNFA (right panel) as well as IL10 and IL5 (not shown) by FACS-ELISA (Miltyni). Values were quantile normalized to reduce batch effects. Following ANOVA, Tukey's multiple comparison test was performed between Mock and each gene target to determine statistical significance. (*) $p < 0.05$; (**) $p < 0.01$; (***) $p < 0.005$; (****) $p < 0.0001$.

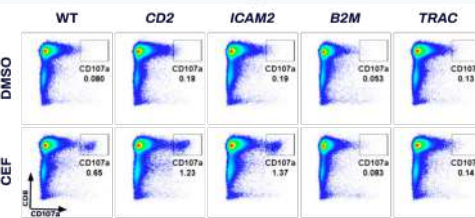

CD8+ T cell antigen-specific cytotoxic activity. Edited CD8+ T cells were thawed and stimulated for 10 days with HLA-I matched, activated B cells loaded with CMV, EBV and IAV peptide mix (CEF). Antigen re-stimulated cells were incubated for 4 hours with CEF peptides in the presence of FITC labeled anti-CD107a antibody. Levels of CD107a expression in gated live, CD3+ cells were compared to unedited control (Mock). As expected, antigen specific responses were completely abrogated in TRAC and B2M edited CD8+ cells, compared to Mock, CD2 or ICAM-2 edited cells. These results are consistent with the canonical roles of B2M and TRAC.